



Luna CPAP & Auto-CPAP Patient Set-Up Manual

1. Button Functions



Mute Button: Silences the alarm and displays iCode® strings




Dial: Starts CPAP pressure and adjusts device settings



Ramp Button: Press it to enable the Ramp feature

2. Setting up the Device

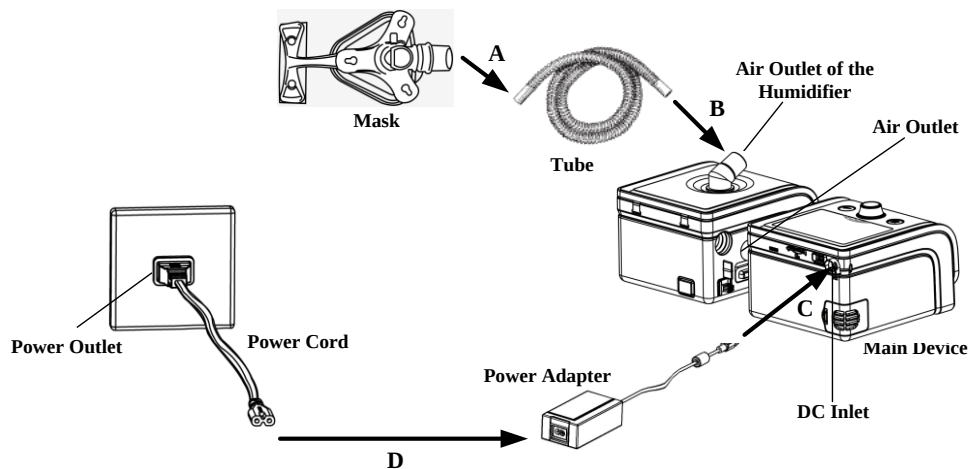

A. Connect one end of the tube to the assembled mask.

B. Connect the other end of the tube securely to the air outlet of the main device or of the humidifier, as shown in the figure above.



C. Connect the power adapter to the DC inlet of the main device.

D. Connect the power cord to the power adapter.

E. Plug the power cord into the power outlet.

3. Adjusting Device Settings

A. After the device is connected to power, its screen displays the Main Interface shown in the left figure below (only applies to Luna Auto-CPAP), or the Main Interface shown in the right figure below (only applies to Luna CPAP).

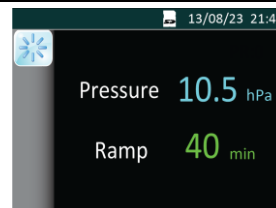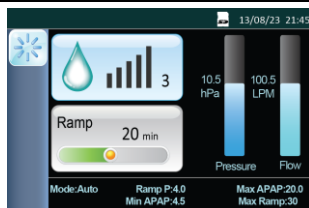


B. From the Main Interface shown above, press and hold the **Ramp Button**  for three seconds. The screen displays the Initial Setup Interface, as shown in the figure below.



C. Press **the Dial**  to access the Setup Interface, as shown in the figure below.



D. Press and turn **the Dial**  to adjust an option.

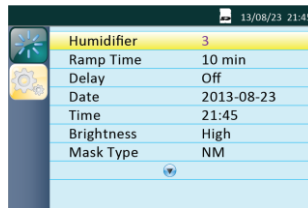

E. Press **the Dial**  to confirm your adjustment to the option.



F. The cursor scrolls through the options as you turn **the Dial** . Refer to steps C-E to adjust other parameters. When the cursor is on the **Home** option, as shown in the figure below, press **the Dial**  to exit the Setup Interface.



4. Starting Treatment

After you finish adjusting device settings and return to the Main Interface, press **the Dial**  to start treatment.

5. Important Tips

Please refer to the user manual for any information on the use of the SpO₂ Kit and SD card.

Before using the device, please read the entire user manual and memorize all information necessary to ensure proper use of the device.